

LOW IRON KIT



Code 7787-01 | Low Range Comparator, 0.05-1.0 ppm

QUANTITY	CONTENTS	CODE
30 mL	*Iron Reagent #1	*4450-G
4.5g	*Iron Reagent #2 Powder	*4451-S
3	Test Tubes, glass, 10mL, w/caps	0843
1	Spoon, 0.15g, plastic	0727
1	Low Range Comparator Viewer	1102
1	Iron Low Range Comparator Bar, 0.05-1.0 ppm	7788-01



*Reagent is a potential health hazard. **READ SDS:** lamotte.com

Emergency information:
Chem-Tel USA 1-800-255-3924
Int'l, call collect, 813-248-0585




To order individual reagents or test kit components, use the specified code number.

Warning! This set contains chemicals that may be harmful if misused. Read cautions on individual containers carefully.
Not to be used by children except under adult supervision.

USE OF THE LOW RANGE COMPARATOR

The Low Range Comparator allows extremely faint colors to be matched to color standards by viewing the reaction down the length of the tube increasing the path length by ten times compared to comparators that view the reaction across the diameter of the tube.

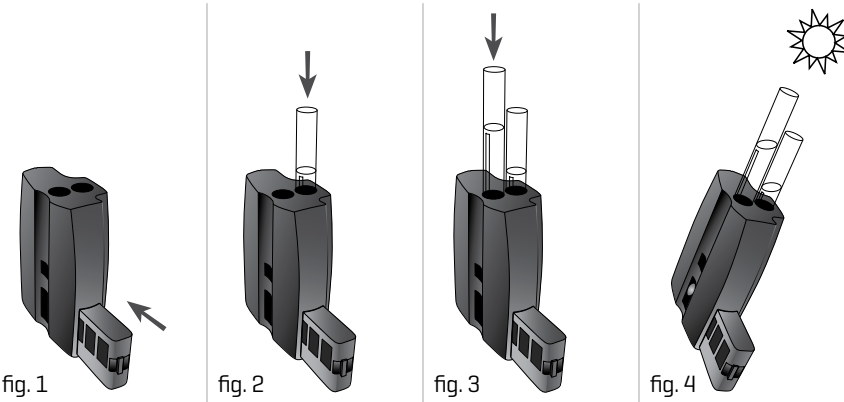

1. Slide the Low Range Comparator Bar into the Low Range Comparator Viewer [fig.1].
2. Fill a test tube to the 10 mL line with untreated sample water. Insert it in the rear hole on the top of the Low Range Comparator [fig.2].
3. Fill a test tube with sample water. Follow the test procedure. Remove the cap. Insert the tube in the front hole on the top of the Low Range Comparator [fig.3].
4. Position the comparator so that light shines down through the test tubes [fig.4]. Tilt the comparator until the color standards and sample are illuminated. Match the color of the reaction to the color standards. Read the result from the Low Range Comparator Bar.

PROCEDURE

1. Slide the Iron Low Range Comparator Bar [7788-01] into the Low Range Comparator Viewer [1102]. Fill a test tube [0843] to the 10mL line with sample water. Place in Low Range Comparator.
2. Add ten drops of *Iron Reagent #1 [4450]. Cap and invert to mix.
3. Use the 0.15g spoon [0727] to add one level spoonful of *Iron Reagent #2 [4451]. Cap and mix until dissolved. Wait 3 minutes.
4. Remove the cap from the test tube. Place the tube in the Low Range Comparator. Match the sample color to a color standard. Record as ppm iron.