

DC1500 UDV ABSORBANCE



UDV Absorbance Colorimeter Kit

Code 3250

CONTENTS

DC1500 UDV Absorbance 568 nm Kit	
■ Kit Contents	
Accessories	3
Analysis	
■ Procedure	3
Calibration	
■ Procedure	5
Set Up	
Setting the Clock	
■ Logging Data	
■ Factory Setup	8
Setting Power Save	
Setting Backlight Time	
Selecting A Language	11
■ Looping Menus	11
Computer Connection	
Output	
■ Computer Connection	12
Battery	
■ Battery/AC Operation	12
Battery Replacement	13
Maintenance	
■ Cleaning	13
■ Returns	13
■ Meter Disposal	14
General Operating Information	
Overview	
■ General Operating Information	14
■ The Keypad	
Sample Holders	
■ The Display and Menus	
■ Tubes and Chambers	16
Sample Dilution Techniques	17
General information	
Packaging and Delivery	
■ General Precautions	
Safety Precautions	17
■ Limits of Liability	
■ Specifications	18
■ CE	18
■ IP67	18
■ Warranty	19
Trouble Shooting Guide	
■ Error Messages	19
Stray Light	
■ Troubleshooting Guide	19

DC1500 UDV ABSORBANCE

■ KIT CONTENTS

QUANTITY	CONTENTS	CODE
4	Square Vials, 10 mm	29653
4	Vial Caps	28687
1	Syringe, 3 mL	1184
1	Water Sample Collecting Bottle	0688
1	1500 Colorimeter for UDV Absorbance, 568 nm	27926-UDV
1	USB Wall Adapter	1721
1	USB Cable	1720-01
1	DC1500-UDV Absorbance 568 nm Colorimeter, Manual	3250-MN
1	DC1500-UDV Absorbance 568 nm Colorimeter, Quick Start Guide	3250-QG

*Reagent is a potential health hazard. **READ SDS:** lamotte.com **Emergency information:**





Chem-Tel USA 1-800-255-3924 | Int'l, call collect, 813-248-0585

ACCESSORIES


DESCRIPTION	CODE
Replacement Chamber	3-0038
COD/UDV Adapter	1724
USB Cable	1720-01
USB Wall Adapter	1721
Car Charger	5-0132

WARNING: Only use the USB Cable (1720-01) that is supplied with the kit. Make no substitution.

ANALYSIS

DDUCEDIIDE

FN	JUEDUKE			
1.	Fill the Water Sample Collecting Bottle (0688) with sample water.			
2.	Press and hold until colorimeter		Main Menu	
	turns on.	Testing Men	и	
		Editing Men	u	
		12:00:00	001/500	(IIIII)
3.	Press ENTER to select Testing Menu .		All Tests	
		103 Absorba	ance 568	
		12:00:00	001/500	(IIIII

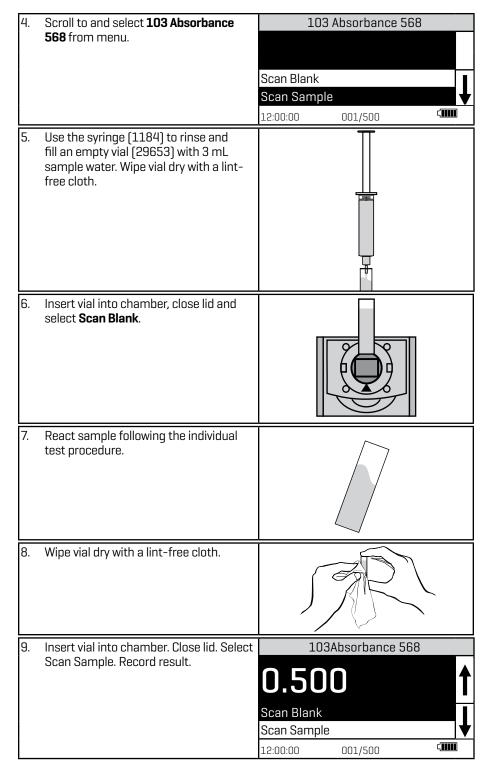


NOTE: The meter will remember the last scanned blank reading. It is not necessary to scan a blank each time the test is performed. To use the previous blank reading, instead of scanning a new one, scroll to Scan Sample and proceed. For the most accurate results, the meter should be blanked before each test and the same tube should be used for the blank and the reacted sample.


NOTE: For the best possible results, carry a reagent blank through the procedure. After scanning the blank, perform the test procedure on clear, colorless, distilled or deioized water. Subtract results of regent blank from all subsequent test results.


CALIBRATION

The calibration should be done with a distilled or deionized water blank and one standard of known concentration. The concentration of the calibration standard should be similar to the expected concentration of the sample that will be tested.

PROCEDURE

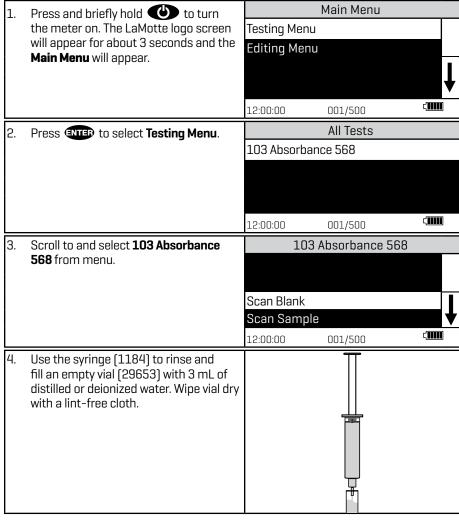

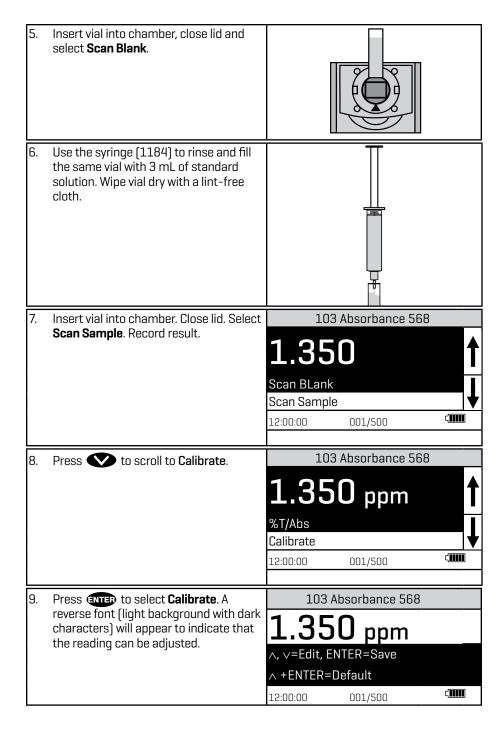


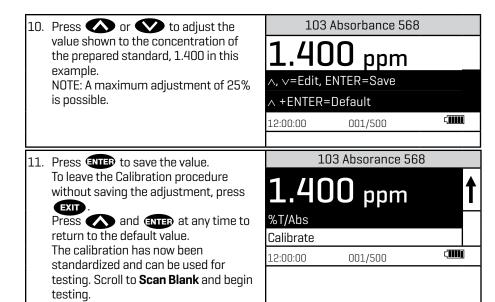

SET UP

■ SETTING THE CLOCK

Setting the clock allows the correct time and date stamp to be stored with each reading in the data logger.

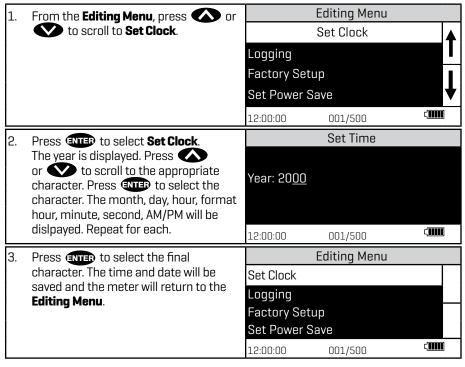

■ LOGGING DATA

The default setting for the data logger is enabled. The meter will log the last 500 data points. The counter in the center bottom of the display will show how many data points have been logged. The display will show 500+ when the data logger has exceeded 500 points and the data points are being overwritten.


Pon	its and the data points are being overwrite	011.		
1. From the Editing Menu , press or			Editing Menu	
	to scroll to Logging .	Set Clock		A
		Logging		
		Factory Se	etup	
		Set Power	Save	▼
		12:00:00	001/500	(IIII)
2.	Press ENTER to select Logging.		Logging	
		Display Tes	t Log	
		Logging Er		
		Logging Di	isabled	
		Erase Log		
<u></u>		12:00:00	001/500	(IIII)
3.	Press or to scroll to		Logging	
	desired function.	Display Tes	st Log	
		Logging En		
		Logging Di	isabled	
		Erase Log		<u> </u>
Ļ		12:00:00	001/500	
4.	Press ENTER. The screen will display	0 : 01 - 1	Editing Menu	
	Storing for about 1 second and return to the Editing Menu .	Set Clock		→
		Logging		I
		Factory Se	etup	
		Set Power	Save	—
		Set Fower	Jave	

■ FACTORY SETUP

The Factory Setup menu is used in manufacturing of the colorimeter. This menu is not for use by the operator in the field.


■ SETTING POWER SAVE

The power saving Auto Shutoff feature will turn the meter off when a button has not been pushed for a set amount of time. The default setting is disabled. To change the setting:

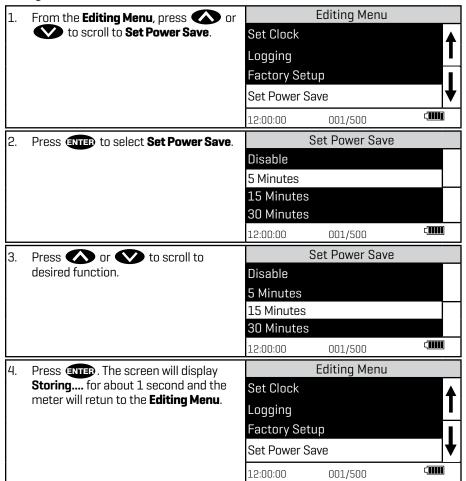


■ SETTING THE BACKLIGHT TIME

The backlight illuminates the display for enhanced viewing. The default setting is 10 seconds. If Button Control is chosen the backlight button on the key pad will act as an on/off switch and the backlight will remain on or off when the meter is being used. When one of the other settings – 10, 20 or 30 seconds – is chosen, the display will be illuminated for the specified amount of time after any button is pressed. NOTE: The backlight feature uses a significant amount of power. The longer the backlight is on, the more frequently the battery will have to be charged if the USB/Wall Adapter is not being used.

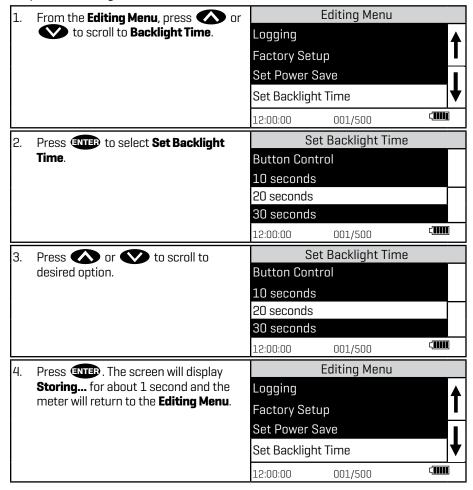

SELECTING A LANGUAGE

There are eight languages available: English, Spanish, French, Portuguese, Italian, Chinese, Japanese, and Turkish.

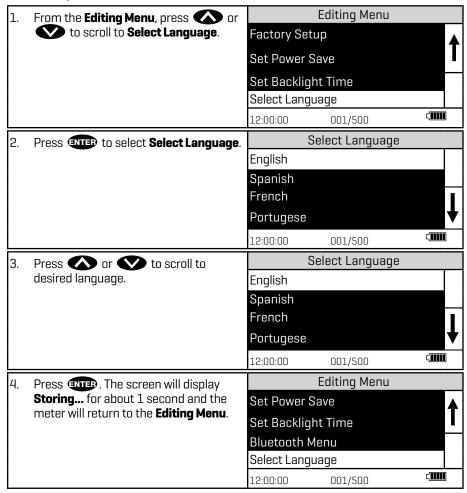

NOTE: If meter unintentionally switches to another language, use the procedure above to reset the meter to the desired language. For example, to reset the meter to English: Turn meter on.

- 1. Press one time. Press
- 2. Press five times. Press ITER.
- 3. Press ENTER.

■ LOOPING MENUS

Long menus, such as the Tsting Menu, incorporate a looping feature which allows the user to quickly reach the last choice in the menu from the first choice. In a looping menu the last choices in the menu are above the first choice and scrolling upward moves through the menu in reverse order. Scrolling downward moves through the menu from first choice to last but the menu starts over following the last choice. So all menu choices can be reached by scrolling in either direction. The diagrams below demonstrate a looping menu.

AND SO ON	AND SO ON	AND SO ON
: : :		: : :
: : :	: : :	
THIRD TO LAST	LAST CHOICE	LAST CHOICE
SECOND TO LAST	AND SO ON	AND SO ON
LAST CHOICE	: : :	: : :
TESTING MENU	TESTING MENU	TESTING MENU
FIRST CHOICE	.	: : : .
SECOND CHOICE	THIRD TO LAST	THIRD TO LAST
THIRD CHOICE	SECOND TO LAST	SECOND TO LAST
ANOTHER	LAST CHOICE	LAST CHOICE
12:00:00 001/500 c	12:00:00 001/500	12:00:00 001/500 c
AND ANOTHER	FIRST CHOICE	FIRST CHOICE
AND SO ON	SECOND CHOICE	SECOND CHOICE
: : :	THIRD CHOICE	THIRD CHOICE
1 : :	ANOTHER	ANOTHER
LAST CHOICE	AND ANOTHER	AND ANOTHER

The feature called **Looping Menu** can be turned on and off in the **Editing Menu**. The default setting is ON.

COMPUTER CONNECTION

USB

■ COMPUTER CONNECTION

USB Type A, USB mini B, Order Cable Code 1720-01.

■ WATERLINK CONNECT 2

The meter may be interfaced with any Windows-based 64-bit computer by using the LaMotte WaterLink Connect 2 program and a USB cable. The program will store test information and results in a local database and allow for exporting this data to a comma separated value [CSV] file. The meter will send the following data: name of test, sample value, sample units, time/date stamp, and meter name. To download WaterLink Connect 2 go to softwarecenter.lamotte.com. Select "WaterLink Connect 2 FREE Download".

BATTERY

■ BATTERY/AC OPERATION

The DC1500 may be operated on battery power using a USB wall adapter or USB computer connection. If using the meter as a bench top unit, use the wall adapter if possible to extend the battery life. The meter will remain on when the USB adapter is used.

To charge the lithium ion battery with the wall adapter, plug the smaller end of the USB cable [USB mini B connector] into the meter and the larger end of the USB cable [USB type A connector] into the wall adapter. Plug the wall adapter into an AC outlet. Reinsert the rubber USB port plug after charging. Failure to insert the USB port plug when the meter is not connected to a computer by USB or actively charging by USB could result in damage to internal components.

To charge the battery from a computer, plug the smaller end of the USB cable [USB mini B connector] into the meter and the larger end of the USB cable [USB Type A connector] into a USB port on a computer. Reinsert the USB port plug after charging. Failure to insert the USB port plug when the meter is not connected to a computer by USB or actively charging by USB could result in damage to internal components.

The battery icon will show no bars and flash when the unit first turns on. Then the indicator will indicate the battery status by showing 0, 1, 2, 3, or 4 bars. It will take 6 hours to fully charge a low battery. The battery icon will flash when the battery is charging. The battery icon will show four bars and stop flashing when it is fully charged. The charging circuit will automatically switch to a float charge when the battery is fully charged. The charger may remain connected. Some computers will NOT supply power to their USB ports during standby operation. The wall adapter will charge the unit continuously. Storing the meter above ambient room temperature will decrease the battery charge more quickly than storage at room temperature. If the meter does not turn on, it means that the battery is at a very low charge. Charging the battery with the wall adapter in this state may take up to 10 hours. At low temperatures, approaching 0 °C, the battery will charge more slowly. It will not charge at all below 0 °C.

The battery icon will show no bars and continuously flash if the battery is getting low but the unit will still operate normally. A "Low Battery" message on the status bar of the display will replace the time when the battery voltage is too low for proper operation and accuracy may be degraded. A "Shutdown Low Batt" message on the display will appear for a few seconds before the power is switched off when the battery is too low to operate the unit. When the battery icon simultaneously flashes bars 1 and 2 followed by bars 3 and 4, it is an indication that the battery is damaged and technical support should be contacted.

To extend the life of the battery:

- Shut down the unit with the power switch when not taking measurements or use the power save option to have the unit automatically turn off after 5 minutes.
- Store the unit at a moderate temperature.
- Fully charge the battery before storing the unit for extended periods of time.
- Fully charge the battery at least once per year. Failure to do so may result in a permanently drained battery.
- Limit backlight use. The unit consumes three times the normal power when the backlight is on. Set the backlight time option to 10 seconds or select "Button Control" and keep the backlight off.

■ BATTERY REPLACEMENT

The lithium ion battery used in this unit should last for many years with normal use. When it no longer powers the unit long enough to meet testing requirements it should replaced. Contact LaMotte Company by phone [1-800-344-3100] or email [tech@ lamotte.com] for more information.

MAINTENANCE

CLEANING

Clean the exterior housing with a damp, lint-free cloth. Do not allow water to enter the light chamber or any other parts of the meter. To clean the light chamber and optics area, point a can of compressed air into the light chamber and blow the pressurized air into the light chamber. Use a cotton swab dampened with Windex* window cleaner to gently swab the interior of the chamber. Do not use alcohol; it will leave a thin residue over the optics when dry.

■ RETURNS

Should it be necessary to return the meter, pack the meter carefully in a suitable

container with adequate packing material. A return authorization number must be obtained from LaMotte Company by calling 800-344-3100 (US only) or 410-778-3100, faxing 410-778-6394, or emailing tech@lamotte.com. Often a problem can be resolved over the phone or by email. If a return of the meter is necessary, attach a letter with the return authorization number, meter serial number, a brief description of problem and contact information including phone and FAX numbers to the shipping carton.

METER DISPOSAL

Waste Electrical and Electronic Equipment (WEEE). Natural resources were used in the production of this equipment. This equipment may contain materials that are hazardous to health and the environment. To avoid harm to the environment and natural resources, the use of appropriate take-back systems is

recommended. The crossed out wheeled bin symbol on the meter encourages the use of these systems when disposing of this equipment.

Take-back systems will allow the materials to be reused or recycled in a way that will not harm the environment. For more information on approved collection, reuse, and recycling systems contact local or regional waste administration or recycling services.

GENERAL OPERATING INFORMATION

OVERVIEW

The meter is a portable, microprocessor controlled, direct reading colorimeter. It has a graphical liquid crystal display and 6 button keypad. These allow the user to select options from the menu driven software, to directly read test results or to review stored results of previous tests in the data logger. The menus can be displayed in eight languages.

The LaMotte tests are precalibrated for LaMotte reagent systems. The colorimeter displays the result of these tests directly in units of concentration.

The optics feature a colored LED. The LED has a corresponding silicon photoiode with an integrated interference filter. The interference filter selects a narrow band of light from the corresponding LED for the colorimetric measurements.

A USB wall adapter, USB computer connection or lithium battery powers the colorimeter.

■ GENERAL OPERATING INFORMATION


The operation of the colorimeter is controlled by the menu driven software and user interface. A menu is a list of choices. This allows a selection of various tasks for the colorieter to perform, such as scan blank and scan sample. The keypad is used to make menu selections that are viewed on the display.

■ THE KEYPAD

The keypad has 6 buttons which are used to perform specific tasks.

	This button will scroll up through a list of menu selections.		
ENTER	The button is used to select choices in a menu viewed on the display.		
	This button controls the backlight on the display.	ENTER	-0-
	This button will scroll down through a list of menu selections.	EXIT	
EXIT	This button exits to the previous menu.		
	This button turns the meter on or off.		

SAMPLE HOLDERS

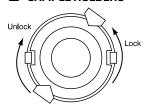


The sample chamber is designed for 25 mm round tubes. An adapter to hold 10 mm square UDV cuvettes is included. Position the COD/UDV Adapter [Code 1724] so that the notches in the adapter fit around the posts on the chamber. Turn the adapter counterclockwise until the arrows are at the top and bottom of the chamber and the adapter is locked into place. Turn the adapter clockwise to unlock the adapter and remove it from the chamber.

■ THE DISPLAY AND MENUS

The display allows menu selections to be viewed and selected. These selections instruct the colorimeter to perform specific tasks. The menus are viewed in the display using two general formats that are followed from one menu to the next. Each menu is a list of choices or selections.


The display has a header line at the top and a footer line at the bottom. The header displays the title of the current menu. The footer line displays the time and the date, the data logger status and the battery status. The menu selection window is in the middle of the display between the header and the footer.

The menu selection window displays information in two general formats. In the first format only menu selections are displayed. Up to 4 lines of menu selections may be displayed. If more selections are available they can be viewed by pressing the arrow buttons to scroll the other menu selections into the menu selection window. Think of the menu selections as a vertical list in the display that moves up or down each time an arrow button is pressed. Some menus in the colorimeter are looping menus. The top and bottom menu choices are connected in a loop. Scrolling down past the bottom of the menu will lead to the top of the menu. Scrolling up past the top of the menu will lead to the bottom of the menu.

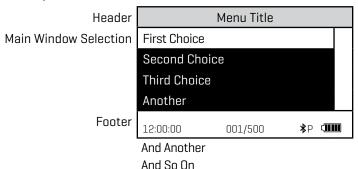

A light bar will indicate the menu choice. As the menu is scrolled through, the light bar will highlight different menu choices. Pressing the the button will select the menu choice that is indicated by the light bar.

In the second format the menu choice window takes advantage of the graphical capabilities of the display. Large format graphic information, such as test results or error messages or the LaMotte logo is displayed. The top two lines of the display are used to display information in a large, easy to read format. The menus work in the same way as previously described but two lines of the menu are visible at the bottom of the display.

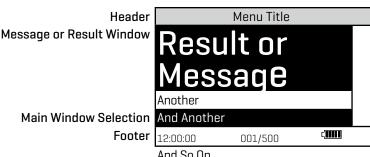

And So On Last Choice

As described previously, the total button allows an exit or escape from the current menu and a return to the previous menu. This allows a rapid exit from an inner menu to the main menu by repeatedly pushing the turn the colorimeter off.

The display may show the following messages:

(IIIII)	Battery Status
†	More choices are available and can be viewed by scrolling up and/or down through the display.
Header	Identifies the current menu and information on units and reagent systems if applicable.
Footer	In the data logging mode the number of the data point is displayed and the total number of data points in the memory will be shown. The footer also shows current time and battery status

■ TUBES AND CHAMBERS

The handling of the tubes is of utmost importance. Tubes must be clean and free from lint, fingerprints, dried spills and significant scratches, especially the central zone between the bottom and the sample line.

Scratches, fingerprints and water droplets on the tube can cause stray light interference leading to inaccurate results. Tubes that have been scratched in the light zone through excessive use should be discarded and replaced with new ones.

Tubes should always be washed on the inside and outside with mild detergent prior to use to remove dirt or fingerprints. The tubes should be allowed to air-dry in an inverted position to prevent dust from entering the tubes. Dry tubes should be stored with the caps on to prevent contamination.

After a tube has been filled and capped, it should be held by the cap and the outside surface should be wiped with a clean, lint-free absorbent cloth until it is dry and smudge-free. Handling the tube only by the cap will avoid problems from fingerprints. Always set the clean tube aside on a clean surface that will not contaminate the tube. It is imperative that the tubes and light chamber be clean and dry. The outside of the tubes should be dried with a clean, lint-free cloth or disposable wipe before they are placed in the meter chamber.

Tubes should be emptied and cleaned as soon as possible after reading a sample to prevent deposition of particulates on the inside of the tubes.

Variability in the geometry of the glassware and technique is the predominate cause of variability in results. Slight variations in wall thickness and the diameter of the tubes

may lead to slight variations in the test results. To eliminate this error the tubes should be placed in the chamber with the same orientation each time.

Chambers which have been scratched through excessive use should be discarded and replaced with a new one.

■ SAMPLE DILUTION TECHNIQUES

If a test result is out of the range of the meter, it must be diluted. The test should then be repeated on the diluted sample. The following table gives quick reference guidelines for dilutions of various proportions.

Amount of Sample	Deionized Water to Bring Final Volume to 10 mL	Multiplication Factor
10 mL	0 mL	1
5 mL	5 mL	2
2.5 mL	7.5 mL	4
1 mL	9 mL	10
0.5 mL	9.5 mL	20

All dilutions are based on a final volume of 10 mL, so several dilutions will require small volumes of the water sample. Graduated pipets should be used for all dilutions. If volumetric glassware is not available, dilutions can be made with the colorimeter tube. Fill the tube to the 10 mL line with the sample and then transfer it to another container. Add 10 mL volumes of deionized water to the container and mix. Transfer 10 mL of the diluted sample to the colorimeter tube and follow the test procedure. Repeat the dilution and testing procedures until the result falls within the range of the calibration. Multiply the test result by the dilution factor. For example, if 10 mL of the sample water is diluted with three 10 mL volumes of deionized water, the dilution factor is four. The test result of the diluted sample should be multiplied by four.

GENERAL INFORMATION

■ PACKAGING AND DELIVERY

Experienced packaging personnel at LaMotte Company assure adequate protection against normal hazards encountered in transportation of shipments.

After the product leaves LaMotte Company, all responsibility for safe delivery is assured by the transportation company. Damage claims must be filed immediately with the transportation company to receive compensation for damaged goods.

■ GENERAL PRECAUTIONS

READ THE INSTRUCTION MANUAL BEFORE ATTEMPTING TO SET UP OR OPERATE THE

METER. Failure to do so could result in personal injury or damage to the meter. The meter should not be used or stored in a wet or corrosive environment. Care should be taken to prevent water from wet tubes from entering the meter chamber.

NEVER PUT WET TUBES IN THE METER.

■ SAFETY PRECAUTIONS

Read the labels on all LaMotte reagent containers prior to use. Some containers include precautionary notices and first aid information. Certain reagents are considered hazardous substances and are designated with a * in the instruction manual. Safety Data Sheets (SDS) can be found at www.lamotte.com. Read the SDS before using these reagents. Additional emergency information for all LaMotte reagents is available 24 hours a day from the National Poison Control Center 1-800-222-1222 or by contacting the 24 hour emergency line for ChemTel 1-800-255-3924 (USA, Canada, Puerto Rico); locations outside the North American Continent 813-248-0585 (call collect). Be prepared to supply the name and four-digit LaMotte code number found on the

container label or at the top of the SDS or in the contents list of the procedure. LaMotte reagents are registered with a computerized poison control information system available to all local poison control centers.

Keep equipment and reagent chemicals out of the reach of young children.
Ensure that the protection provided by this equipment is not impaired. Do not intall or use this equipment in a manner that is not indicated in this manual.

■ LIMITS OF LIABILITY

Under no circumstances shall LaMotte Company be liable for loss of life, property, profits, or other damages incurred through the use or misuse of its products.

■ SPECIFICATIONS & RANGES

INSTRUMENT TYPE: Colorimeter			
Readout	160 x 100 backlit LCD, 20 x 6 line graphical display		
Wavelengths	568 nm		
Wavelength Accuracy	±2% FS		
Readable Resolution	0.01 Absorbance Units		
Wavelength Bandwidth	10 nm typical		
Photometric Range	-2 to +2 AU		
Photometric Precision	± 0.001 AU at 1.0 AU		
Photometric Accuracy	±0.005 AU at 1.0 AU		
Sample Chamber	10 mm vials		
Light Sources	1 LED		
Detectors	1 silicon photodiode		
Modes	Absorbance		
Languages	English, Spanish, French, Portuguese, Italian, Chinese, Japanese, Turkish		
Temperature	Operation: 0-50 °C; Storage: -40-60 °C		
Operation Humidity Range	0-90 % RH, non-condensing		
USB Port	Mini B		
Power Requirements	USB wall adapter, USB computer connection or lithium ion rechargeable battery		
Battery	Charge Life: Approximately 380 tests with backlight on to 1000 tests with backlight off. Battery Life: Approximately 500 charges.		
Electrical Rating	Rated voltage (5V), Rated power of input current (1.0A) at mini-USB input port		
Data Logger	500 test results		
Waterproof	IP67 with USB port plug in place		
Dimensions (LxWxH)	3.5 x 7.5 x 2.5 inches, 8.84 x 19.05 x 6.35 cm		
Weight	13 oz, 362 g (meter only)		

■ CE COMPLIANCE

The colorimeter has been independently tested and has earned the European CE Mark of compliance for electromagnetic compatibility and safety. To view certificates of compliance, go to the LaMotte website at www.lamotte.com.

■ IP67 Certification

The meets IP67 standards for protection against dust and immersion only when the USB port plug is in place. Documentation is available at www.lamotte.com.

■ WARRANTY

LaMotte Company warrants this instrument to be free of defects in parts and workmanship for 2 years from the date of shipment. Keep the proof of purchase for warranty verification. If it should become necessary to return the instrument during or the warranty period, contact our Technical Service Department at 1-800-344-3100 or 1-410-778-3100, ext. 3 or softwaresupport@lamotte.com for a return authorization number or visit www.lamotte.com for troubleshooting help. The sender is responsible for shipping charges, freight, insurance, and proper packaging to prevent damage in transit. This warranty does not apply to defects resulting from action of the user such as misuse, improper wiring, operation outside of specification, improper maintenance or repair, or unauthorized modification. LaMotte Company specifically disclaims any implied warranties or merchantability or fitness for a specific purpose and will not be liable for any direct, indirect, incidental, or consequential damages. LaMotte Company's total liability is limited to repair or replacement of the product with a new or refurbished meter as determined by LaMotte Company. The warranty set forth above is inclusive and no other warranty, whether written or oral, is expressed or implied.

TROUBLESHOOTING

■ ERROR MESSAGES


Over Range

If the message OVERRANGE is displayed when scanning a sample, the sample may be over range or under range. If the sample is over range, the sample should be diluted and tested again (see Sample Dilution Techniques and Volumetric Measurements, page 17).

If overrange is displayed, scroll to **Calibrate** and press extent to view the result. The concentration is only an **approximation**. Press and continue testing.

Note: If is not pressed the meter will calibrate on the over range result.

Note: Over range results will also be displayed in the test log.

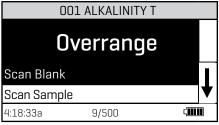

STRAY LIGHT

The colorimeter should have no problems with stray light. Make sure that the sample compartment lid is always fully closed.

■ TROUBLESHOOTING GUIDE

PROBLEM	REASON	SOLUTION
Flashing	Low battery. Readings are reliable.	Charge battery or use USB wall/computer adapter.
"Low Battery"	Battery voltage is very low. Readings are not reliable.	Charge battery or use USB wall/computer adapter.
"Shut Down Low Batt" Shut Down	Battery is too low to operate the unit.	Charge battery or use USB wall/computer adapter.
"Overrange"	Sample is outside of acceptable range.	Dilute sample and test again.
Unusually large negative or positive readings when performing calibration	Incorrect standards used to calibrate meter.	Use fresh 0.0 standard in clean tube. Reset meter to factory default settings. Recalibrate meter.



802 Washington Avenue | Chestertown, Maryland | 21620 USA 800.344.3100 | +1 410.778.3100 (Outside USA) | Fax 410.778.6394 | lamotte.com